



# ASKINA® PAD | ASKINA® PAD S

ABSORBING NON-WOVEN COMPRESS

B. Braun Wound Care. Right. From the start.

# Askina® Pad | Askina® Pad S

## Absorbing non-woven compress

Askina® Pad is a sterile, non adherent, absorptive compress to cover all types of exuding wounds.

Askina® Pad S is a sterile, non adherent, absorptive slit compress for the dressing of draining and feeding tubes. It has a central hole suitable for the use with tubes with diameters up to 11mm and a slit, to wrap around tubing.

- Absorbency capacity: min 5g/25cm<sup>2</sup>
- Spread of blood: min 1.5g/25cm<sup>2</sup>



## INDICATIONS

### Askina® Pad:


- Absorptive dressing for traumatic, dermatologic or surgical wounds

### Askina® Pad S:

- Absorptive dressing for draining sites, around, catheters, feeding tubes, suprapubic catheters, etc.

## THREE LAYER SYSTEM

1. Water repellent backing material
2. Natural cotton fleece with a high absorption capacity
3. Non adherent wound contact layer (polyester net)

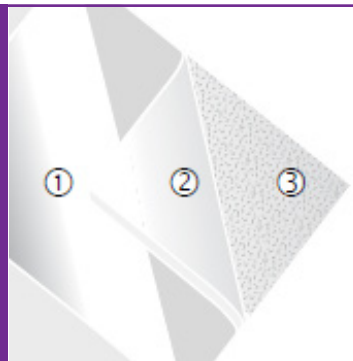

## B. Braun Ordering Information

### Askina® Pad | Askina® Pad S

Askina®	Size	Unit/Box	Code
Pad	5 x 5 cm	100	9024018
Pad	10 x 10 cm	100	9024034
Pad	20 x 10 cm	100	9024042
Pad S	5 x 5 cm	30	9024050
Pad S	7.5 x 7.5 cm	30	9024069