




INTERVENTIONAL  
VASCULAR  
DIAGNOSTICS  
AND THERAPY

## SURECAN<sup>®</sup> SAFETY II

POWER-INJECTABLE SAFETY NEEDLES FOR ACCESS PORTS

# Surecan® Safety II needles for Access Ports

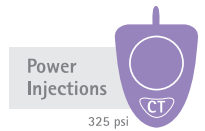

## Surecan® Safety II is a non-coring safety needle with a small size and unique design.

The Surecan® Safety II is easy to handle with maximal safety to minimize needlestick injuries. This provides high comfort for clinicians and patients, in hospitals as well as for home care treatment.

Indicated for up to seven days for:

- Infusion of chemotherapy
- Total parenteral nutrition (TPN)
- Antibiotic therapy
- Transfusion
- Pain relief and drainage
- Blood sampling
- Power injections of contrast media up to 325 psi

### Ergonomic low profile wings

- Provide secure handling during puncture and removal
- Permit small and comfortable dressing
- Needle gauge indicated on wings

### Closed-cell pad

- Increases patient comfort

### Special needle tip

- Aids needle retention in the septum during power injections

### Unique bevel

- Non-coring bevel designed to minimize septum damage

### Coloured clamp

- Possible interruption of infusion at any time
- Easy to identify gauge of the needle
- Contains maximum psi information for high pressure injections

### Extension tubing

- Transparent anti-kink tubing
- Power injectable

Extensive choice of references:  
19 G to 22 G  
12 mm to 38 mm



Surecan® Safety II, compatible with all B. Braun Access Port Systems

## Surecan® Safety II complies with the requirements of the European Council Directive 2010/32/EU from May 10th, 2010: "Prevention from sharp injuries"

- "Health and safety of workers is paramount and is closely linked to the health of patients."\*
- Surecan Safety II is a "medical (device) ... incorporating safety-engineered protection mechanisms."\*


\* European council directive 2010/32/EU from May 10th, 2010

### Green point and "click" confirms safety mechanism

- Easy confirmation of correct deployment
- Double affirmation through visual and audible features




### Transparent base

- Easier visualisation of the insertion point
- Holes aid skin aeration

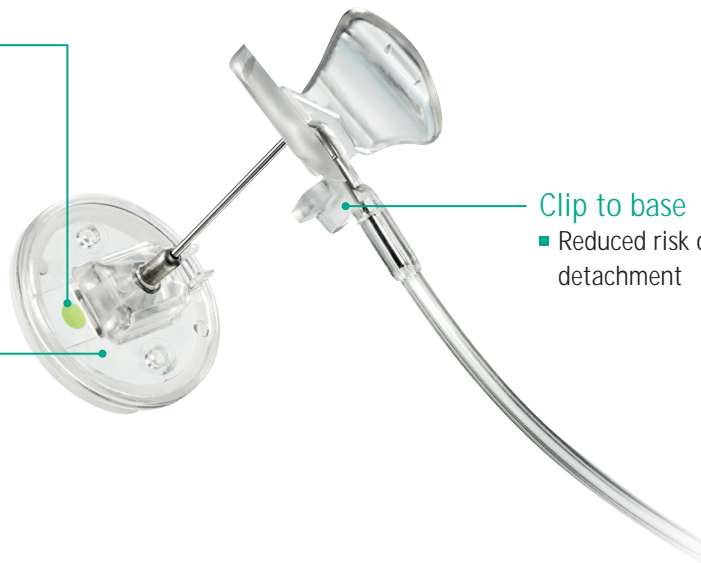

### Clip to base

- Reduced risk of accidental detachment

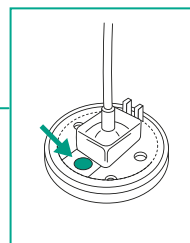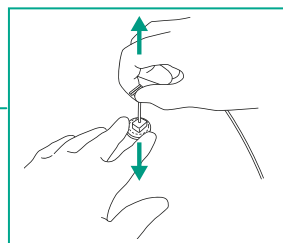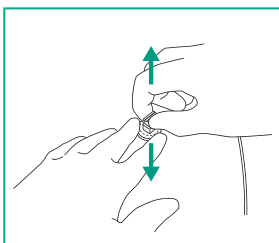

### Safe removal


- Stabilise the base of the needle
- Firmly pull the wings up until you hear a „click“ and see the green point

For more information on the risks and prevention of needle-stick injuries, please visit [www.safeinfusiontherapy.com](http://www.safeinfusiontherapy.com) and refer to our brochure no. 6069095.

## Caresite® as pre-connected needle-free valve


The Surecan® Safety II is also available with the pre-connected needle-free valve Caresite®. The split septum valve technology combined with the ergonomic design enables clear and easy access.

- Protects against microbiological contamination
- Prevents catheter occlusion through positive displacement (valve design)
- Ease of use

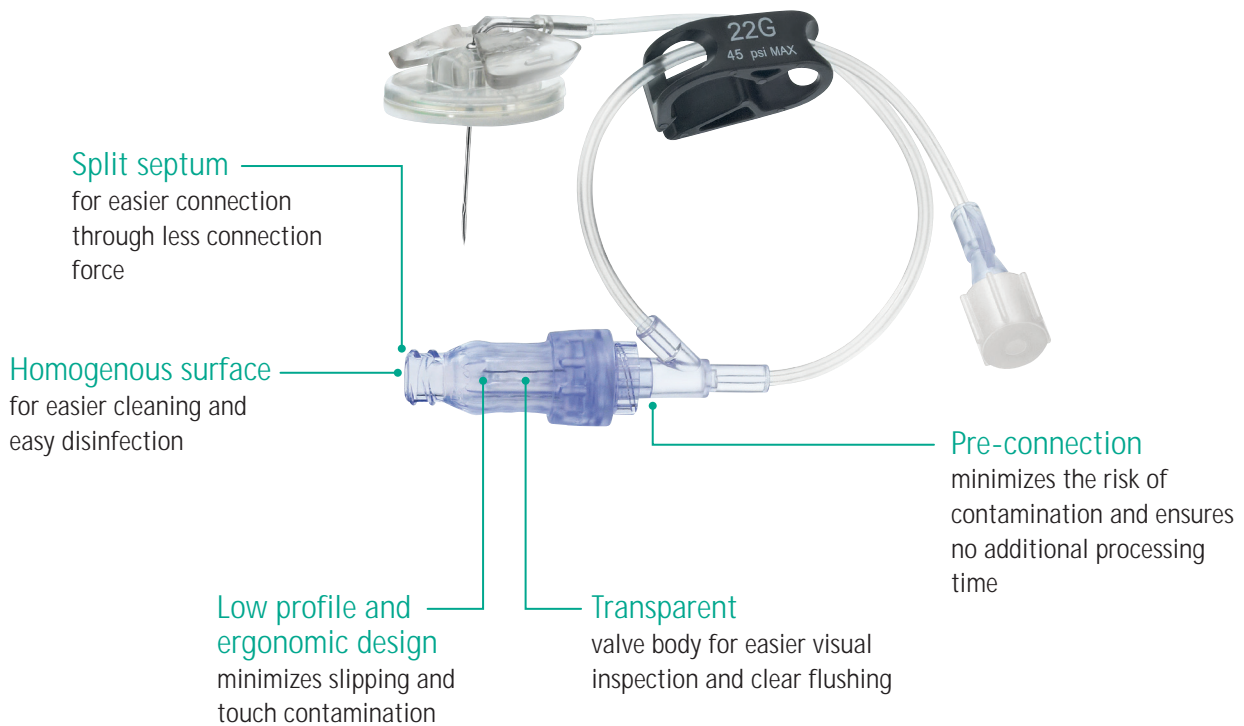


### ATTENTION

High pressure (HP) injections have to be performed through the Caresite® valve

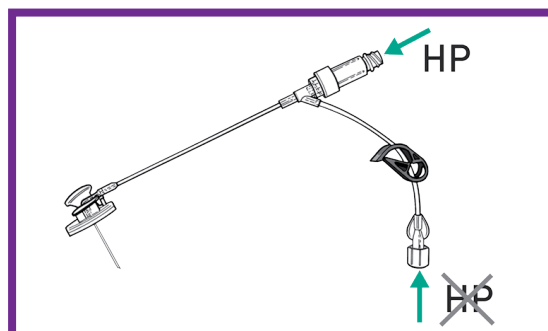


## Surecan® Safety II safety needle

### Surecan® Safety II

	Size	Pri <sup>^</sup> Z_X Volume/mL	Cannula Length (mm)	Unit	Reference								
 <table border="1" data-bbox="135 705 406 884"> <tr> <td>POWER</td> <td>MR</td> </tr> <tr> <td>INJECTABLE</td> <td>CONDITIONAL</td> </tr> <tr> <td>DEHP</td> <td>LATEX</td> </tr> <tr> <td>FREE</td> <td>FREE</td> </tr> </table>	POWER	MR	INJECTABLE	CONDITIONAL	DEHP	LATEX	FREE	FREE	G19	0.30	12	20 per box	04447042
	POWER	MR											
	INJECTABLE	CONDITIONAL											
	DEHP	LATEX											
	FREE	FREE											
	15	04447000											
	20	04447001											
	25	04447002											
	32	04447003											
	38	04447004											
	G20	0.25	12	20 per box	04447043								
			15		04447005								
			20		04447006								
			25		04447007								
			32		04447008								
			38		04447009								
	G22	0.20	12	20 per box	04447044								
			15		04447010								
20			04447011										
25			04447012										
32			04447013										

### Surecan® Safety II with Y-site and Caresite®

	Size	Pri <sup>^</sup> Z_X Volume/mL	Cannula Length (mm)	Unit	Reference								
 <table border="1" data-bbox="135 1467 406 1646"> <tr> <td>POWER</td> <td>MR</td> </tr> <tr> <td>INJECTABLE</td> <td>CONDITIONAL</td> </tr> <tr> <td>DEHP</td> <td>LATEX</td> </tr> <tr> <td>FREE</td> <td>FREE</td> </tr> </table>	POWER	MR	INJECTABLE	CONDITIONAL	DEHP	LATEX	FREE	FREE	G19	0.65	12	20 per box	04447057
	POWER	MR											
	INJECTABLE	CONDITIONAL											
	DEHP	LATEX											
	FREE	FREE											
	15	04447045											
	20	04447046											
	25	04447047											
	32	04447048											
	38	04447049											
	G20	0.55	12	20 per box	04447058								
			15		04447050								
			20		04447051								
			25		04447052								
			32		04447053								
			G22		0.45	12	20 per box	04447059					
	15	04447054											
	20	04447055											
25	04447056												

B. Braun Australia Pty Ltd | Bella Vista NSW | Tel.1800 251 705 | info.au@bbraun.com | www.bbraun.com.au  
B. Braun New Zealand | Customer Care 0800 657 775 | Fax (09) 526 0782 | www.bbraun.co.nz

BANZ AE E0397 Rev. A 08/18 | ©2018 B. Braun Australia